


SDS-PAGE Standard Protein Marker (Middle Range)

ChemWhat Technical Data Sheet (TDS)

Catalog Number:	801-01
Quantity:	20 lanes/vial
Contents:	Mixed six kinds of lyophilized (freeze-dried) protein powder (40 µg each). By Laemmli
	method, after SDS-PAGE electrophoresis and Coomassie Blue R-250 staining, it presents
	six uniform bands.
Storage:	-20 °C
Recommend Gel	SDS-PAGE 12 %
Concentration:	
Reconstitution:	We recommend that this vial be briefly centrifuged prior to opening to bring the contents to
	the bottom. Reconstitute in 200 μl reduced 1 \times SDS loading buffer and divide to 20 aliquots.
	Store at \leq -20 °C. Melt at room temperature and heat at boiling water bath for 3-5 minutes
	before using.

Design Sketch:

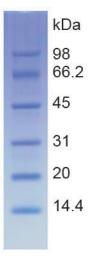

CHEMWHAT LIMITED

https://www.chemwhat.com

Email: contact@chemwhat.com